

Set de Verificación de Calibración para HI 706

HI706-11



imagen no
disponible

Descripción

El set de verificación de Calibración Fosforo Rango Alto Checker® HC - **HI 706-11** proporciona una solución simple para validar el Fosforo Rango Alto Checker® HC - **HI 706** en el rango de 0,00 a 15,00 ppm (mg/L). Este conjunto de estándares de alta calidad se fabrica en nuestras instalaciones de última generación y se suministra con un certificado de análisis. El Certificado de análisis proporciona el número de lote, los valores de referencia y la fecha de vencimiento para la trazabilidad al certificar el Checker® HC.

Especificaciones

No Especifica

Accesorios

No Especifica

Cómo pedir

No Especifica

Ventajas

El **HI 706-11** es un conjunto de estándares de calibración de alta calidad que permite a los usuarios validar el fósforo alto rango Checker® HC - **HI 706**. El **HI 706-11** incluye una cubeta sellada para cero ppm de fósforo y una cubeta teñida a una muestra de fósforo equivalente a 7,5 +/- 0,5 ppm a 25 °C.

El HI 706-11 se suministra con Certificado de Análisis

- Código de producto
- Numero de lote
- Fecha de caducidad



CERTIFICATE OF ANALYSIS

Product Name: CALIBRATION CHECKING SET FOR HI 711

Product Code: HI 711-11

Lot number: SC0372/14

Expiration date: October 2017

Standard cuvette	Lot number	Standard Value @ 25°C (ppm)	Passed
A, HI 711-11	2038	1.00 ± 0.05	✓
B, HI 711-11	2039	1.00 ± 0.05	✓

Hanna Instruments certifies these standards meet the stated tolerance limits as indicated. These should be used for quality control to validate instrument function on the specified Hanna instrument only. Store standards in the storage box upright at room temperature. Do not use standard cuvettes.


Instructions for use:

1. Thoroughly wipe standard cuvette with HI731318 lint free cloth and only handle using top of cuvette. The cuvette wall must be dry and free of smudges.
2. At the C.1 prompt place A, HI 711-11 standard cuvette into the holder of the Checker and press operations button.
3. At the C.2 prompt place B, HI 711-11 standard cuvette into the holder and press the operations button.

Acceptable results must be between 0.95 and 1.05 ppm @ 25°C.

File number: CERT711-11_SC0372/14

QA manager: Bogdan Munteanu

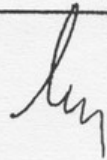

Video

No Especifica